

Gymna ShockMaster 300

Radial shockwave therapy - anytime, anywhere

Product data sheet



User friendly

- Easy to transport thanks to the custom made carry bag



- Silent air compressor, **no maintenance, no emptying of the condensation bottle**
- **Unique automatic condensation**
- Compressor switches on only during treatment
- Rapid software update via USB port



Always up to date

Shockwave therapy is evolving constantly. Thanks to treatment software updates, you stay abreast of new developments:

- **Unique! patient Guided Therapy System:** the intelligent Shockmaster patient Guided Therapy System will automatically adapt treatment level, based on your patient's condition
- **Total therapy plan:** every application is divided in different steps: pGTS, tendon treatment, trigger point treatment, muscle treatment with D-Actor, muscle treatment with V-Actor
- Treatment parameters are adapted in function of the latest information
- New applications are added, as well as relevant treatment protocols




Intuitive user interface

- **Easy application search:** submenu's based on the affected tissue



- **Unique:** marked treatment area and referred pain with trigger point

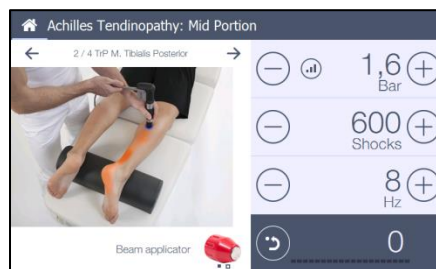

- **Safety restriction:** preventing the use of an inappropriate handpiece
- **Pop-up** after 1.25 million shocks, when revision is needed
- **Contra indications:** easy check list before each treatment

NEW Handpiece Short

- Ergonomic design: lightweight and well balanced
- Integrated soft grip materials
- Extended energy range: starting level at 0,3bar_{eff}
- Broader application: more patients, more applications
- Compatible with new Fascia and Spine applicators
- Backwards Compatibility: the Handpiece Short may be used on all ShockMaster devices provided that the mandatory software update is performed



Fascia applicator set



Spine applicator set

Gymna Club

Free subscription on the Gymna club:

Extra information about shockwave therapy with the Gymna ShockMaster: literature, technical support, software updates for free, extra marketing material

Technical specifications

Characteristics

Dimensions (H x D x W cm)	18 x 38 x 42
Weight (kg)	11
Compressor	Integrated air compressor, silent, vibration free No maintenance with the automatic condensation
Mobile	With the carrier bag
Connections	1
Touch screen	7inch Full color 180° view

User interface

Plug and play	Choose your own settings
List of applications	43 applications
Total therapy plan	Every application is divided in different steps
Pre-programmed settings	> 90 programs
Body part	Choose applications via a body area
Marked treatment area	Treatment area marked on the pictures during treatment Trigger point referred pain application
Quick start	Tendon applications Muscle applications Bone applications Connective tissues Neurological applications
Pre-set counter	Supported via illuminated bar
Software update	Via USB port

Performance

Frequency range	0,5 - 17 Hz (in steps of 0,1 Hz)
Max frequency with V-Actor	31 Hz
Pressure range	0,3 - 4 bar _{eff} (in steps of 0,1 bar _{eff})

Handpiece and applicators

On/off function	Electronic No air cooling
Weight (kg)	0,48 kg
Ergonomics	Integrated soft rubber protection
V-Actor vibration therapy	Available
Revision kit (standard)	4 revisions Do it yourself possibility
Applicators (standard)	Beam applicator D-Actor 20 mm applicator Deep Impact applicator
Applicators (optional)	15mm applicator 6mm applicator D-Actor 35 mm applicator Focus lens applicator

